

# ApoE ToolSet™ for LightCycler™ (ApoE E2/E3/E4, codons 112 Cys/Arg & 158 Cys/Arg) Application on LightCycler 480 II

Lyophilized ToolSet for PCR using the LightCycler™ Instrument. Licensed by Roche Diagnostics GmbH

**Order#: ApoE - 16**

1 ToolSet for 16 reactions

Store at 4°C, protected from light.  
 Exposure to light may especially damage  
 the Oligotool™ tube (vial with red cap).

For use with **LightCycler FastStart DNA Master HybProbe**, 10 x conc. (Roche Cat.No.: 03003248001)

**1. ToolSet contents**

Vial	Label	Content	Quantity
			<b>ApoE - 16</b>
<b>1, Red cap</b>	<b>OligoTool</b>	- lyophilized oligos for PCR - contains mutation detection and anchor probe, primers	For 16 tests  Dissolved: 50 µL
<b>2, Green cap</b>	<b>Control</b>	- lyophilized ApoE E2/E4 HET DNA	Dissolved: 20 µL
<b>3, Blue cap</b>	<b>Solvent</b>	- to dissolve OligoTool / Control	1000 µL of Solvent
<b>4, Yellow Cap</b>	<b>DMSO</b> <b>Caution Irritant !</b> <b>May be absorbed</b> <b>through skin !</b>	- PCR additive for GC rich regions	100 µL of DMSO  Note : Fluid above 18 °C Solid below 18 °C

Additional equipment and reagents required but not supplied :  
 LightCycler FastStart DNA Master HybProbe, 10 x conc. (Roche Cat.No.: 03003248001), including 25mM MgCl<sub>2</sub>;  
 LightCycler instrument, LightCycler capillaries, DNA extraction materials

**2. Introduction**

**2.1. Product overview**

**ToolSet description** This ToolSet is specifically designed for multiplex genotyping of the codon 112 Cys / Arg and codon 158 Cys / Arg polymorphisms in the Apolipoprotein E (ApoE) gene, encoding the **E2/E3/E4 isoforms of Apolipoprotein E**, by LightCycler PCR with Melting Curve Analysis. Primers and fluorescent probes have been optimized for specific amplification of a 462 bp segment containing the relevant site and for optimal genotype discrimination.

**Control material** E2 / E4 Heterozygous Control DNA, lyophilized  
 (= Cys112Arg / Cys158Arg compound heterozygote).

**Storage of ToolSet and Solutions** Store at +4°C when lyophilized, protected from light.  
 The unopened lyophilized ToolSet is stable at +4°C for 12 months from date of manufacture if protected from light. When dissolved store at +4°C for a maximum of 4 weeks, or at -20°C for longer periods (months), protected from light. Avoid freezing and thawing > 3 times.

### 3. Preparation for LightCycler PCR

**Toolset preparation** **Dissolve** the content of the **OligoTool** tube (Red Cap) with **50 µl of Solvent**.  
**Dissolve** the content of the **Control** tube (Green Cap) with **20 µl of Solvent**.

1. Before opening tubes, centrifuge them quickly.
2. Add Solvent into OligoTool tube and Control tube as above.
3. Recap tubes, vortex gently.
4. Before opening tubes, centrifuge them quickly.
5. Proceed to Reaction Mix preparation.

**Primers ?** You don't have to add primers.

**Probes ?** You don't have to add probes.

**Reaction Mix Preparation** For 1 (One) reaction, prepare the Reaction Mix as shown in the following table :

Reagent	µL
OligoTool ApoE -16, dissolved	2.8
Solvent ApoE -16	8.0
MgCl <sub>2</sub> 25 mM	1.2 (final 2.5 mM)
DMSO (yellow cap)	2.0
FastStart DNA Master HybProbe, 10x	2.0
Total Reaction Mix	16.0
+ Your DNA or Control ApoE E2/E4 HET -16	4.0
Grand Total	20.0

Use FastStart DNA Master HybProbe 10x and MgCl<sub>2</sub> 25 mM from Roche LightCycler FastStart DNA Master HybProbe, 10 x conc. (Roche Cat.No.: 03003248001, including 25mM MgCl<sub>2</sub>).

For multiple reactions, multiply the indicated volumes appropriately.

**Positive Control** Always run a positive control with the samples. Use the dissolved heterozygous Control ApoE E2/E4 HET - 16 DNA (Green Cap).

**Negative control** Always run a negative control with the samples. To prepare a negative control, replace the template DNA with Solvent (Blue Cap).

**Extraction of genomic DNA** You can use different Kits for DNA isolation, either with a manual method or with an automated system. The elution buffers should be salt-free. Example : Roche High Pure PCR Template Preparation Kit (Cat.No. 1 796 828)

**Application** The **ApoE ToolSet™** for LightCycler™ allows detection of the **codon 112 Cys/Arg and codon 158 Cys/Arg variants** in the ApoE gene, encoding the **E2/E3/E4 isoforms of Apolipoprotein E**.

Note : This ToolSet was developed for use in life science research only.

#### 4. LightCycler 480 Settings and Experimental Protocol

For use with LC 480 Program Version 1.5 series.

Detection : Dynamic Red **640** (498-640 nm) and Dynamic Red **705** (498-660 nm)  
with color compensation (obligatory !)

##### Fast Start Enzyme Activation and DNA Denaturation

Cycle Program Data	Value
Cycles	1
Analysis Mode	None
Temperature Targets	<b>Segment 1</b>
Target Temperature (°C)	<b>95</b>
Incubation time (s) <b>Hold</b>	<b>600</b>
Temperature Transition Rate (°C/s) <b>Ramp Rate</b>	<b>4.4</b>
<b>Acquisitions (per °C)</b>	---
Secondary Target Temperature (°C)	0
Step Size (°C)	0
Step Delay (Cycles)	0
Acquisition Mode	None

##### Amplification

Cycle Program Data	Value		
Cycles	<b>45</b>		
Analysis Mode	Quantification		
Temperature Targets	<b>Segment 1</b>	<b>Segment 2</b>	<b>Segment 3</b>
Target Temperature (°C)	<b>95</b>	<b>60</b>	<b>72</b>
Incubation time (s) <b>Hold</b>	<b>5</b>	<b>10</b>	<b>15</b>
Temperature Transition Rate (°/s) <b>Ramp Rate</b>	<b>4.4</b>	<b>2.2</b>	<b>4.4</b>
<b>Acquisitions (per °C)</b>	---	---	---
Secondary Target Temperature (°C)	0	0	0
Step Size (°C)	0	0	0
Step Delay (Cycles)	0	0	0
Acquisition Mode	None	<b>Single</b>	None

##### Melting Curve Analysis

Cycle Program Data	Value		
Cycles	<b>1</b>		
Analysis Mode	Melting Curves		
Temperature Targets	<b>Segment 1</b>	<b>Segment 2</b>	<b>Segment 3</b>
Target Temperature (°C)	<b>95</b>	<b>40</b>	<b>85</b>
Incubation time (s) <b>Hold</b>	<b>30</b>	<b>120</b>	---
Temperature Transition Rate (°/s) <b>Ramp Rate</b>	<b>4.4</b>	<b>1.5</b>	<b>0.1</b>
<b>Acquisitions (per °C)</b>	---	---	<b>3</b>
Secondary Target Temperature (°C)	0	0	0
Step Size (°C)	0	0	0
Step Delay (Cycles)	0	0	0
Acquisition Mode	None	None	<b>Continuous</b>


##### Cooling

Cycle Program Data	Value
Cycles	<b>1</b>
Analysis Mode	None
Temperature Targets	<b>Segment 1</b>
Target Temperature (°C)	<b>40</b>
Incubation time (s) <b>Hold</b>	<b>30</b>
Temperature Transition Rate (°/s) <b>Ramp Rate</b>	<b>1.5</b>
Secondary Target Temperature (°C)	0
Step Size (°C)	0
Step Delay (Cycles)	0
Acquisition Mode	None

## 5. Typical results

### Introduction

Use the Melting Curve program to genotype the human genomic DNA research samples. The melting peaks allow discrimination between the possible genotypes of the codon 112 and 158 polymorphisms in the ApoE gene. Figure 1 shows a typical result obtained with the ApoE ToolSet™ for LightCycler™ :

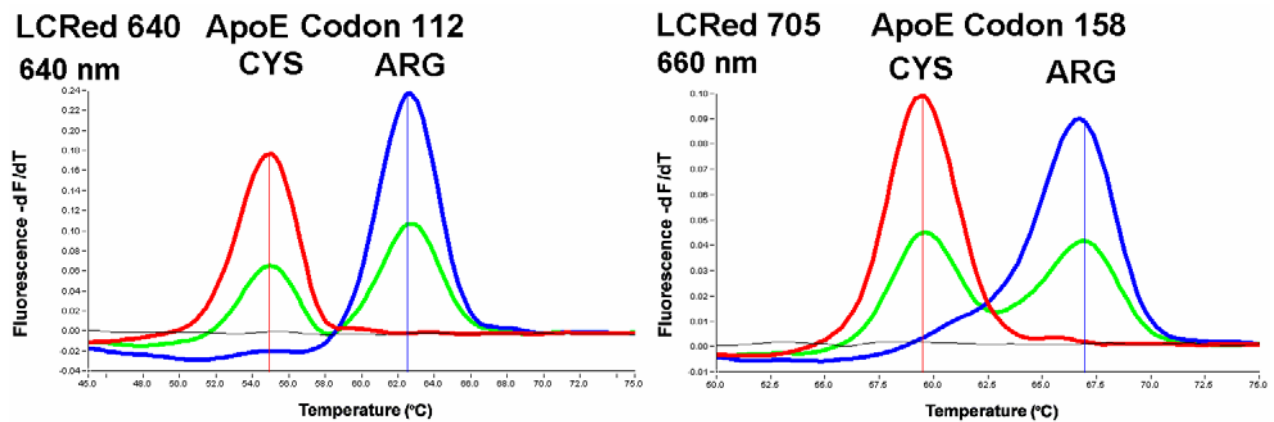

**Figure 1 : Melting curve analysis of genotypes of codons 112 and 158 in the ApoE gene.**

**BLUE :** Homozygotes **Arg112Arg** (left) and **Arg158Arg** (right)  
**GREEN :** Compound Heterozygote **Cys112Arg / Cys158Arg** = ApoE E2 / E4 HET Control contained in the ToolSet.  
**RED :** Homozygotes **Cys112Cys** (left) and **Cys158Cys** (right)  
**BLACK :** No DNA Control.

Conditions : LC 480 program version 1.5, With Color compensation,

**For conversion of codon 112 / 158 melting curves into ApoE genotypes, use the following Table :**

**Table 1 : Conversion of Codon 112 / Codon 158 melting curve results into ApoE genotype**  
**YES : peak observed at indicated Tm      no : No peak observed at indicated Tm**

Genotype	Codon 112 (LCRed 640, 640nm)		Codon 158 (LCRed 705, 660 nm)		Genotype
	Tm 54 - 56 C	Tm 62 - 64 C	Tm 59 - 61 C	Tm 66 - 68 C	
E2/E2 (Cys/Cys)	YES	no	YES	no	E2/E2 (Cys/Cys)
E3/E3 (Cys/Arg)	YES	no	no	YES	E3/E3 (Cys/Arg)
E4/E4 (Arg/Arg)	no	YES	no	YES	E4/E4 (Arg/Arg)
E2/E3	YES	no	YES	YES	E2/E3
E2/E4	YES	YES	YES	YES	E2/E4
E3/E4	YES	YES	no	YES	E3/E4

**Note :** The values for the respective melting temperatures may vary for +/- 2.5 °C between different experiments. The Delta T between the melting peaks for different genotypes may vary +/- 1.0 °C. The ApoE ToolSet™ has been developed for and validated with the LightCycler™ and its original accessory materials and reagents. Performance of the ToolSet with other instruments, accessories and reagents has not been validated by ratiogen.

### 7. Notices to Purchaser : Licenses and Trademarks, Prohibition of Resale

ratiogen's ToolSets™ for LightCycler™ are produced and marketed under license from Roche Diagnostics GmbH. The purchase of this product does not convey any right for its use in clinical diagnostics. No license to use the PCR process is conveyed expressly or by implication by purchase of this product. The LightCycler technology is licensed from Idaho Technology Inc., Idaho Falls, ID, USA. LightCycler and High Pure are trademarks of members of the Roche group. Resale of ratiogen's products is expressly prohibited.

### How to contact ratiogen

E-mail [info@ratiogen.com](mailto:info@ratiogen.com)

Internet <http://www.ratiogen.com>